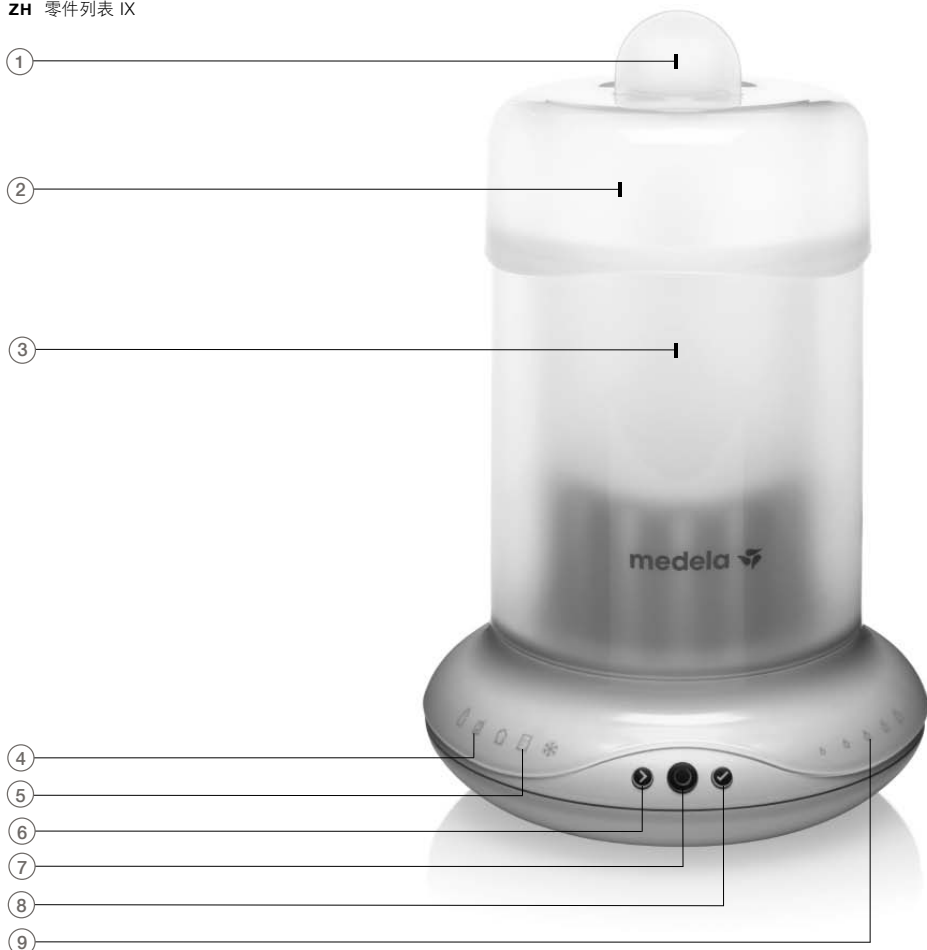
*No. 1 choice
of hospitals
& mothers*

Bottle Warmer

EN Instructions for use
DE Gebrauchsanweisung
FR Mode d'emploi
IT Istruzioni per l'uso
NL Gebruiksaanwijzing
ES Manual de instrucciones
PT Manual de instruções
EL Οδηγίες χρήσης
HE הוראות שימוש
SV Bruksanvisning
DA Brugsanvisning
FI Käyttöohjeet
NO Bruksanvisning
PL Instrukcja obsługi
RU Инструкция по применению
HU Használati utasítás
CS Návod na použití
TR Kullanma talimatları
JA 使用説明書
KO 사용 설명서
TW 使用說明書
ZH 使用说明书

EN Parts List IX
DE Teileliste IX
FR Liste des pièces IX
IT Elenco componenti IX
NL Onderdelenlijst IX
ES Lista de piezas IX
PT Lista de peças IX
EL Κατάλογος εξαρτημάτων IX
HE IX רשימת חלקים
SV Lista över delar IX
DA Reservedelsliste IX
FI Osaluettelo IX
NO Reservedelsliste IX
PL Wykaz części zamiennych IX
RU Список деталей IX
HU Alkatrészek listája IX
CS Seznam dílů IX
TR Parça Listesi IX
JA 部品リスト IX
KO 부품 명칭 IX
TW 配件清單 IX
ZH 零件列表 IX

EN	Instructions for use	5
DE	Gebrauchsanweisung	14
FR	Mode d'emploi	23
IT	Istruzioni per l'uso	32
NL	Gebruiksaanwijzing	41
ES	Manual de instrucciones	50
PT	Manual de instruções	59
EL	Οδηγίες χρήσης	68
HE	הוראות שימוש	85
SV	Bruksanvisning	86
DA	Brugsanvisning	95
FI	Käyttöohjeet	104
NO	Bruksanvisning	113
PL	Instrukcja obsługi	122

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Breastfeeding Solution Circle

Dear Parents,

Congratulations on your decision to purchase a Medela product. For over 40 years, Medela has helped to enhance babies' health through the life-giving benefits of breast milk and supported mothers during their breastfeeding experience – making it as safe and as easy as possible. Medela strives to understand nature and over the years has become one of the world's leading experts in breastfeeding.

Breastfeeding is one of the most beautiful and natural experiences you can share with your baby. It's a unique time when you feel intensely close and intimate. Medela provides products for your various needs during breastfeeding. The Medela Bottle Warmer which is easy to transport and clean helps you to safely warm your precious breast milk. Medela wishes you and your baby all the best!

Product Overview

 On/Off	 Select			 Enter	
			 ml/ oz		
 Breastmilk			 Baby Food		
 Room-Temperature ~20°C	 Fridge-Temperature ~7°C		 Freeze-Temperature ~-18°C		
	 30 ml/1.0 oz	 80 ml/2.7 oz	 150 ml/5.0 oz	 200 ml/6.7 oz	 250 ml/8.5 oz
	 < 150 g	 > 150 g	not selectable	not selectable	not selectable

I Bottle Warmer

The new Medela Bottle Warmer brings liquids or food to a constant 34°C – be they frozen, out of the fridge or at room temperature. Helps preserve nutrients and vitamins.

II Indications for Use

- The Medela Bottle Warmer
- heats milk bottles and pre-cooked baby food to the ideal feeding temperature (around 34°C +/- 3°C) and stabilises the temperature.
 - is intended for home use only.

III Important Safeguards

When using electrical devices, basic safety precautions should always be followed.

Read all instructions carefully before use.

Product

- Use the product only for its intended use as described in this instruction manual.
- Only use the power cord that comes with the Bottle Warmer.
- Keep appliance dry and do not place in direct sunlight.
- Never immerse the Bottle Warmer in water.
- Only use the Medela components supplied with the unit.
- Check components for wear or damage.

Electrical

- Make sure the voltage of the device is compatible with the power source.
- Plug the power cord into the Bottle Warmer first and then into wall socket.
- Always unplug electrical devices immediately after use.
- Separation from the mains is only assured through the disconnection of the power cord from the wall socket outlet.
- Keep the power cord away from heated surfaces.
- Never operate an electrical device if it has a damaged power cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.
- Do not reach for any electrical device if it has fallen into water. Unplug immediately.
- If the supply cord is damaged, it must be replaced by a new cord set immediately to avoid a hazard.

Using / Feeding

- This appliance is intended for use by responsible persons, keep out of the reach of children.
- This appliance is not intended for use by persons with reduced physical, sensory or mental capabilities, or lack of experience or knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure they do not play with the appliance.
- There are no user-serviceable parts inside the Bottle Warmer. Do not attempt to repair the Bottle Warmer yourself.
- Always use on a flat, stable and heat-resistant surface.
- Whilst pressing the buttons, please hold the Bottle Warmer base to prevent from slipping.
- Always ensure that the correct volume of water has been added before switching on the Bottle Warmer.
- When in use, the Bottle Warmer contains hot water. Always use and handle the equipment with care.
- Breastmilk/Babyfood should not be heated for too long.

Danger of scalding!

- Take care when opening the lid after a completed cycle as some residual steam may escape.
- Use a cloth when removing bottles or jars, as they may be hot.
- Use the bottle lift to warm small bottles or baby food jars.
- For hygiene and health reasons, once heated, milk and food should not be reheated and should be discarded if not used immediately.
- The final temperature of the food and the heating time could vary according to the volume of food being heated, the initial (room, fridge or freezer) temperature, the consistency of the food, as well as of the type and material of the feeding bottle or jar.
- Always check milk or food temperature before feeding the baby!
Milk: Shake the bottle or storage container to equalise the milk temperature and check the temperature by pouring some drops of liquid on the back of your hand.
Food: Stir the food thoroughly and test it yourself before feeding the baby.

III Important Safeguards

Cleaning

- Failure to follow de-scaling instructions may cause irreparable damage to the product.
- Do not touch the hotplate in the Bottle Warmer. Danger of scalding!
- For further advice and recommendations on baby feeding, please consult your health care professional.
- Always let the Bottle Warmer cool down completely before starting a new warming cycle.
- Waiting time for next use = 15 min.
- If the Bottle Warmer is switched on and not running a cycle, it will automatically turn off after 1 minute.

SAVE THIS INSTRUCTION MANUAL

IV Operation

Warming up feeding bottles & baby food jars

1

Place the appliance on a flat stable surface. Plug the power cord into Bottle Warmer and into wall socket as well. Don't switch on until water and bottles/food jar are in the appliance.

2

Fill the Bottle Warmer with the appropriate amount of water.

Breastmilk				Baby Food	
Milk	Water	Milk	Water	Food	Water
30 ml	120 ml	150 ml	240 ml	< 150 g	120 ml
80 ml	120 ml	200 ml	240 ml	> 150 g	240 ml
		250 ml	240 ml		

Note: Measure with a Medela-Bottle the appropriate amount of water. Using cool distilled or filtered water may reduce limescale build-up. For other amounts of breastmilk, we recommend to choose the next lower cycle for warming the milk, e.g. for 120 ml breastmilk choose:

Milk	Water
80 ml	120 ml

3

Place filled bottle or food jar, in the bottle lift (3) and both into the warmer filled with water.

4

Close the lid (2) of the Bottle Warmer.

5

Press -Button
-Icon blinks. Start your selection process.

6

Select Breast Milk or Baby Food (4)
 To warm up breast milk press Enter-Button to confirm (-Icon turns on)
 To warm up baby food press Select-Button (-Icon blinks)
 Press Enter-Button to confirm (-Icon turns on)

After you have finished your selection, -Icon blinks.

7

Select starting temperature - Room / Fridge / Freezer (5)
 To choose room temperature press Enter-Button to confirm (-Icon turns on)
 To choose fridge temperature press Select-Button
 Press Enter-Button to confirm (-Icon turns on)
 To choose freezer temperature press Select-Button and after
 Enter-Button to confirm (-Icon turns on)

After you have finished your selection, -Icon (30ml / 1.0 oz) blinks.

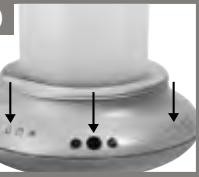
8

Select the required volume (9)
 By using Select-Button , you can switch between different
 volumes (see chart on next page)

To confirm the volume press Enter-Button

9

After you have finished your selection, all selected icons turn on,
-Icon turns orange and heating process has been started.

IV Operation

Breastmilk – Timetable

		 30 ml/1.0 oz	 80 ml/2.7 oz	 150 ml/5.0 oz	 200 ml/6.7 oz	 250 ml/8.5 oz
		3 min	4 min	6 min	7 min	8 min
		5 min	6 min	8 min	10 min	12 min
		18 min	21 min	24 min	not selectable	not selectable

Baby Food – Timetable

		 < 150 g	 > 150 g			
		6 min	7 min	not selectable	not selectable	not selectable
		9 min	10 min	not selectable	not selectable	not selectable
		not selectable	not selectable	not selectable	not selectable	not selectable

Note: Heating up to 34 °C (+/- 3 °C). In order to avoid scalding – do not open the lid while the warmer is in operation.

10

⏻-Button blinks green when the warming cycle is completed (Temperature of 34 °C +/- 3 °C is achieved).

The Bottle Warmer will beep five times (slow). Press the Enter-Button ⏻ to stop.

11

Take the breast milk bottle or food jar with the bottle lift (3) out from the warmer.

(!) If you do not take out the Milk Bottle or Food Jar the temperature will automatically constantly remain at 34 °C (+/- 3 °C) for 20 min.

Always check milk or food temperature before feeding the baby!

12

⏻-Button turns constant green.

After 20 min, Bottle Warmer will beep five times (fast) and the appliance will turn off automatically.

13

Unplug from the mains and empty residual water. If warming an additional bottle or food jar, continue from step 1. Between two warming cycles, always allow the Bottle Warmer to cool down completely (15 min).

The Medela Bottle Warmer can be stored inside the Steam Steriliser.

V Cleaning and Hygiene

After each use

- Always unplug the appliance from the mains before cleaning it.
- Do not clean the appliance when hot.
- Wait 15 minutes after use to allow the Bottle Warmer to cool down.
- Ensure that the plug does not come in contact with water.
- Wipe the inside and outside of the Bottle Warmer with a damp (not wet) cloth.

*Note: Use water of drinking quality for any described cleaning procedures
Do not clean with abrasive cleaners
Never place the Bottle Warmer directly under tap water
Never clean the Bottle Warmer in a dishwasher*

Cleaning

After a while in use limescale will build-up around the hotplate of the Bottle Warmer. This should be removed at regular intervals.

Every month (De-scale)

Version 1

- Mix **100 ml/3.3 oz** of water with citric acid* according to manufacturer's instructions and pour solution into empty Bottle Warmer.
- DO NOT place a bottle or jar into the Bottle Warmer when cleaning.
- DO NOT put lid (2) on.
- Select breastmilk, room temperature and the volume of 250 ml. After you have finished your selection, all selected icons turn on and -Icon turns orange. The cleaning process has been started.
- After the cleaning process has ended, leave the solution to cool for 30 minutes before emptying the Bottle Warmer and rinsing with tap water.
- Wipe the inside and outside of the Bottle Warmer with a damp (not wet) cloth only.

Version 2

- Pour 50 ml/1.8 oz of vinegar mixed with 50 ml/1.8 oz of cold water into the Bottle Warmer.
- DO NOT place a bottle or jar into the Bottle Warmer when cleaning.
- DO NOT start the warming process.
- Leave liquid in the unit until limescale has dissolved.
- Empty the unit and rinse thoroughly with clear water.
- Wipe the inside and outside of the Bottle Warmer with a damp (not wet) cloth only.

**Citric acid is available from most drugstores and pharmacies*

Note: Do not use anything other than citric acid or vinegar as this may cause irreparable damage

VI Troubleshooting

Bottle Warmer will not turn on

Check that the plug is in the Bottle Warmer and in the wall socket also. Then press the -Button.

Bottle Warmer is beeping

If the Bottle Warmer beeps five times (slow), warming cycle is complete.

If the Bottle Warmer beeps five times (fast), 20min of constantly warming will be over and the appliance turns off automatically.

VII Warranty

Medela AG ("Medela"), warrants to the original end user ("Customer") that this product will be free from defects in workmanship and materials, under normal use, for two (2) years.

Instructions are strictly respected, in particular in case of extreme and/or continuous applications/use of the product.

The warranty period begins on the proved purchase date (conclusion of purchase contract). Products purchased from non-authorised distributors, e.g. internet market places, auctions or other types of public sale (including but not limited to eBay) are not covered by this warranty. Medela's sole obligation under this express warranty shall be, at Medela's sole option and expense, to replace or repair the Product or part thereof (in case the product is not available, to deliver a similar type of the Product) or to refund the purchase price paid for the Product. Medela warrants any repaired or replaced Product or part for a period of ninety (90) days from delivery, or through the end of the original warranty, whichever is longer. All Products or parts that are replaced become the property of Medela. This express warranty does neither cover consumables, nor third party products.

Obtaining Warranty Service. Customer must contact the retailer where the product was purchased from, the authorised distributor of Medela or Medela's Service Centre indicated by Medela within the applicable warranty period to obtain warranty service authorisation. Dated proof of original purchase from Medela, its authorised distributor or its authorised retailer and a description of the defect will be required. Medela is not responsible for Products or parts received without a warranty service authorisation. Repaired or replacement Products will be delivered to Customer free of charge. The repaired product or part will be delivered as soon as reasonably possible. Medela shall not be responsible for any damages occurring during such delivery. Medela shall, at its sole discretion, decide on the place of performance for work under warranty. Medela shall have the right to charge for additional costs for services under warranty if the site of the Product is other than where the Product was originally delivered.

Warranty Exclusive. Customer's sole remedy for breach of the warranty shall be the express warranty. The foregoing warranty is exclusive and is in lieu of all other warranties, terms or conditions, expressed or implied, either in fact or by operation of law, statutory or otherwise, including warranties, terms or conditions of merchantability, fitness for a particular purpose, satisfactory quality and non-infringement, all of which are expressly disclaimed. Medela shall not be liable if the alleged defect or malfunction was caused by Customer's or any other person's misuse, neglect, improper handling, unauthorised attempts to open, repair or modify the Product, inadequate maintenance, disregard of operating instructions, excessive load or stress, normal wear and tear, incorrect voltage or any other cause beyond the range of its intended use, by accident, fire, or other hazards, or other cause not due or attributable to Medela. This warranty does not cover physical damage to the Product or malfunctions resulting from the use of the Product in conjunction with any sort of ancillary or peripheral equipment and Medela determines that there is no fault with the Product itself.

Limitation of Liability. Medela also excludes any liability, whether based in contract or tort (including negligence), for incidental, consequential, indirect, special, or punitive damages of any kind, or costs of procurement of substitute products by Customer, or for the loss of revenue or profits, loss of business, loss of information or data, or other information of financial loss arising out of or in consequence with the sale, maintenance, use, performance, failure, or interruption of this product, even if Medela or its distributors have been advised of the possibility of such damages, and limits its liability to replacement, repair, or refund of the purchase price paid, at Medela's option. This Limitation of Liability for damages will not be affected if any remedy provided herein shall fail of its essential purpose.

Disclaimer. Should a court of jurisdiction not allow the entire exclusion or limitation of implied warranties or the limitation of incidental or consequential damages for certain products supplied to consumers or the limitation of liability for personal injury, such implied warranties and such liabilities will be limited to the duration of the applicable express warranty.

Governing Law and Place of Jurisdiction. This Limited Warranty shall be governed by the laws of Switzerland, excluding all conflict of laws principles and excluding the United Nations Convention on the International Sale of Goods dated 11 April 1980. The ordinary courts at the Medela registered office in Baar, Switzerland shall be competent. Medela shall, at its sole discretion, also be entitled to take legal action in the competent courts at the Customer's place of business or domicile. With this Limited Warranty, Medela grants Customer specific legal rights which do not restrict any statutory consumer rights.

VIII Disposal

The unit is made of various metal and plastics. Before disposal, the device is to be rendered unusable and it must not be disposed of as unsorted municipal waste in accordance with local regulations. Use your local return and collection system for waste electrical and electronic equipment. Improper disposal may have harmful effects in the environment and on public health.

IX Parts List

- 1 Lid Handle
- 2 Lid
- 3 Bottle Lift
- 4 Breast Milk or Baby Food Display
- 5 Room-, Fridge- or Freeze-Temperature Display
- 6 Select-Button
- 7 On/Off- and Select/Correct-Button
- 8 Enter-Button
- 9 Volume Display

Note: Whenever you press "on/off"-Button, Bottle Warmer will switch off automatically

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